

C-Myc-Tag Rabbit Polyclonal Antibody

Catalog	TDY023C	TDY023F	Tel: 010-82908854
			Free: 400-0620-621
Quantity	50μL	100μL	Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,	N/A	N/A	IgG


Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**

Recommended dilutions: WB: 1:5,000

Optimal dilutions should be determined by the end user.

Specificity: The C-Myc tag antibody can detect C-Myc fusion proteins.

Background: c-Myc-tag antibody is part of the Tag series of antibodies, the best quality in the research. Myc protein is a transcription factor that activates expression of a great number of genes through binding on consensus sequences (Enhancer Box sequences (E-boxes)) and recruiting histone acetyltransferases (HATs). A recent study demonstrated that temporary inhibition of Myc selectively kills mouse lung cancer cells, making it a potential cancer drug target.

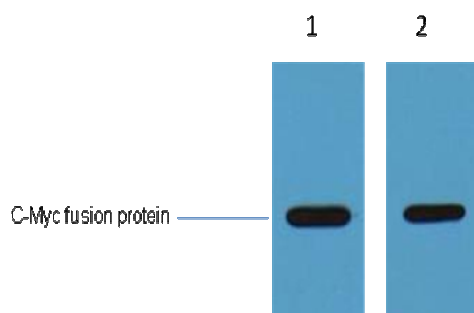

1μg C-Myc fusion protein+ Primary antibody dilution at

1、 1:5,000

2、 1:10,000