



His-Tag Rabbit Polyclonal Antibody

Catalog TDY024C TDY024F

Tel: 010-82908854

Free: 400-0620-621

Quantity 50 μ L 100 μ L

Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,	N/A	N/A	IgG


Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**

Recommended dilutions: WB: 1:5,000

Optimal dilutions should be determined by the end user.

Specificity: The His tag antibody can recognize C-terminal, internal, and N-terminal His-tag fusion proteins.

Background: A **polyhistidine-tag** is an amino acid motif in proteins that consists of at least five histidine (*His*) residues, often at the N- or C-terminus of the protein. Monoclonal antibodies specific to six histidine tags can greatly improve the effectiveness of several different kinds of immunoassays, helping researchers identify, detect, and purify polyhistidine fusion proteins in bacteria, insect cells, and mammalian cells.

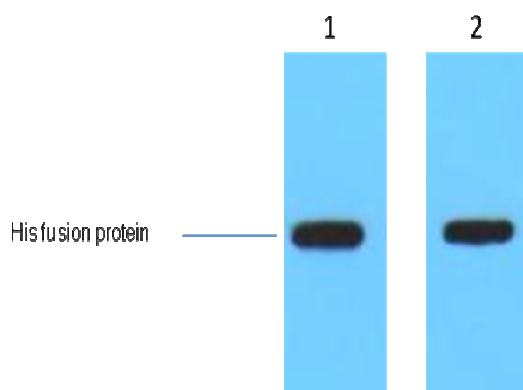

1 μ g His fusion protein+ Primary antibody dilution at

1、 1:5,000

2、 1:10,000