

GFP-Tag Rabbit Polyclonal Antibody

Catalog	TDY026C	TDY026F	Tel: 010-82908854
Quantity	50 μ L	100 μ L	Free: 400-0620-621
			Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,	N/A	N/A	IgG


Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**

Recommended dilutions: WB: 1:5,000

Optimal dilutions should be determined by the end user.

Specificity: The GFP tag antibody can recognize GFP and GFP tag fusion proteins.

Background: The **green fluorescent protein (GFP)** is a protein composed of 238 amino acid residues (26.9kDa) that exhibits bright green fluorescence when exposed to light in the blue to ultraviolet range. In cell and molecular biology, the GFP gene is frequently used as a reporter of expression.

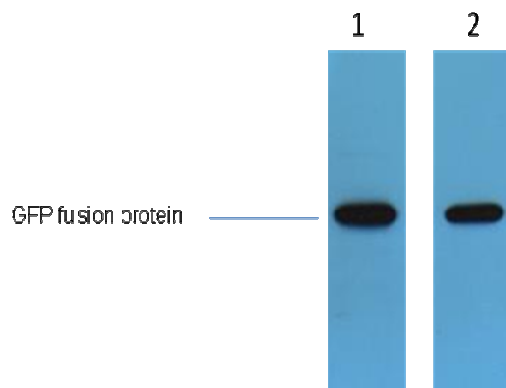

1 μ g GFP fusion protein+ Primary antibody dilution at

1、 1:5,000

2、 1:10,000