

HSV Mouse Monoclonal Antibody(Q72)

Catalog TDY140C TDY140F

Tel: 010-82908854

Quantity 50 μ L 100 μ L

Free: 400-0620-621

Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB	N/A	N/A	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C**. **Do not aliquot the antibody.**


Recommended dilutions: WB: 1:5,000-10,000

Optimal dilutions should be determined by the end user.

Specificity: The HSV tag antibody can recognize HSV tag fusion proteins.

Alternative Names: HSV epitope tag, QPELAPEDPED tag

Background: HSV-tag peptide QPELAPEDPED conjugated to KLH. HSV- Tag antibody can recognize C-terminal, internal, and N-terminal HSV-tagged proteins.

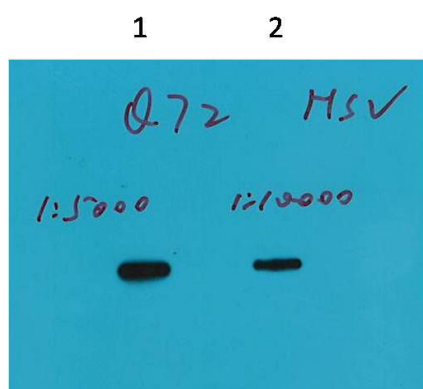

Western blot analysis of Recombinant protein with HSV Mouse mAb(Q72)

diluted at

1) 1:5,000

2) 1:10,000