

SRT-Tag Mouse Monoclonal Antibody(4H9)

Catalog	TDY194C	TDY194F	Tel: 010-82908854
			Free: 400-0620-621
Quantity	50μL	100μL	Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB	N/A	N/A	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C**. **Do not aliquot the antibody.**


Recommended dilutions: WB: 1:5,000

Optimal dilutions should be determined by the end user.

Specificity: The SRT Tag Mouse Monoclonal Antibody can detects SRT Tag recominant protein.

Alternative Names: SRT Tag antibody

Background: The SRT-tag is a peptide epitope comprised of 10 amino acids(TFIGAIATDT) derived from the crystalline surface layer protein of Rickettsia typhi1.

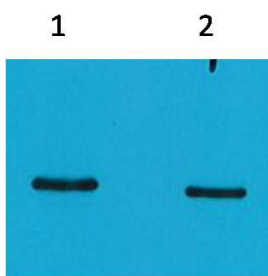

Western blot analysis of SRT-Recombinant protein with TDY194 diluted at

1) 1:5,000

2) 1:10,000