

Rabbit IgG(Fc specific, Heavy Chain) Mouse Monoclonal Antibody(1E1)

Catalog	TDY303C	TDY303F	Tel: 010-82908854
Quantity	50μL	100μL	Free: 400-0620-621
			Web: www.tdybio.com

For research use only.


Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB, IHC	Rabbit	~55KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**

Recommended dilutions: WB: 1:5,000-20,000 IHC: 1:200-500

Optimal dilutions should be determined by the end user.

Specificity: The antibody reacts with the Fc portion of rabbit IgG heavy chain but not with the Fab portion of rabbit immunoglobulins. The antibody may cross-react with immunoglobulins from other species.

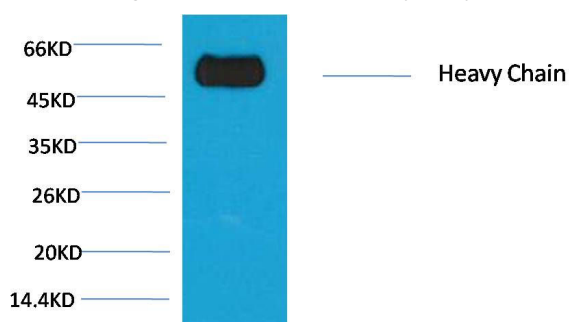

Western blot analysis of Reduced Rabbit IgG, H+L with TDY303 diluted at 1:10,000.