

PI3 Kinase P85 β Rabbit Polyclonal Antibody(F221)

Catalog TDY396C TDY396F

Tel: 010-82908854

Web: www.tdybio.com

Quantity 50 μ L 100 μ L

Entrez-Gene ID: 5296, Swiss-Prot Acc. O00459

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB, IHC	H	82KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C**. **Do not aliquot the antibody.**



Recommended dilutions: WB: 1:1,000-2,000 IHC: 1:50-200

Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous PI3 Kinase P85 β protein.

Alternative Names: P85 antibody, p85 beta antibody, PIK3R2 antibody, Phosphoinositide 3 kinase regulatory subunit 2 antibody

Background: PI3 Kinase P85 β mediating the association of the p110 catalytic unit to the plasma membrane.

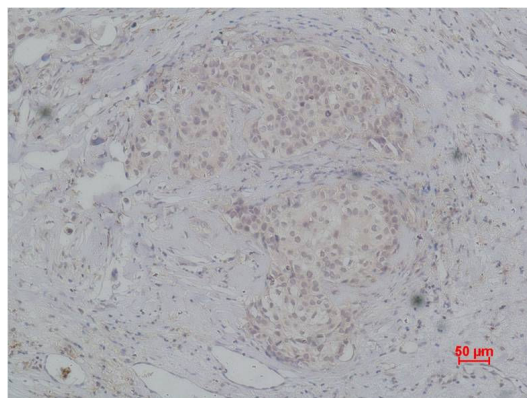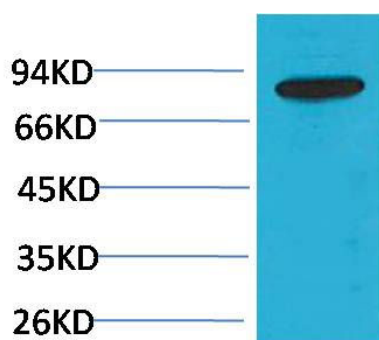

Western blot analysis of HepG2 with PI3 Kinase P85 β Rabbit pAb TDY396 diluted at 1:2,000.

Immunohistochemical analysis of paraffin-embedded human breast carcinoma using PI3 Kinase P85 β (TDY396) Rabbit pAb diluted at 1:200.